

Hand-held retinoscopes for objective refraction.

The BETA 200 Retinoscope is used for exact measurement of refractive error. Both myopia and hyperopia and astigmatism can be measured.

Two versions are available:

The Spot Retinoscope with a round light spot and the more widely-used Streak Retinoscope with a streak image which simplifies the examination.

HEINE ParaStop® and polarisation filter.

The BETA 200 Retinoscope with ParaStop features the latest multi-coated optics for an exceptionally bright fundus reflex and easier detection of the neutralization point. ParaStop was developed by HEINE for the easy, precise selection of a parallel illumination beam. ParaStop simplifies and speeds up the precise detection of cylinder axis. ParaStop also simplifies the verification of the cylinder correction after refraction.

HEINE BETA 200® Retinoscope with HEINE ParaStop®

HEINE ParaStop® for total precision in selecting a parallel beam



- ∴ **ParaStop.** Precise, easy selection of a parallel beam.
- ∴ **Ergonomic shape.** Protects the examiner's orbita from stray light.
- ∴ **XHL Xenon Halogen Technology.** Bright, white light. Bright fundus reflex, easy recognition of neutralization point.
- ∴ **Streak or spot Retinoscope.** Simply by changing the bulb.
- ∴ **Metal controls.** Long-lasting.
- ∴ **Single control for vergence and rotation.** Comfortable operation.
- ∴ **Dustproof.** Maintenance-free.
- ∴ **Integrated polarisation filter.** Eliminates stray light and internal reflexes for a brighter fundus reflex.
- ∴ **Orange filter (optional).** Reduces patient dazzle without affecting the fundus reflex.
- ∴ **Detachable brow rest.** Increased comfort and control during examination.
- ∴ **Holder for fixation cards (optional).** For dynamic retinoscopy.

BETA 200 Streak Retinoscope	2.5V	3.5V
with XHL Xenon Halogen bulb, without handle	C-001.15.353	C-002.15.353
XHL Xenon Halogen spare bulb	X-001.88.087	X-002.88.089
XHL Xenon Halogen spare bulb to convert to spot retinoscope	X-001.88.088	X-002.88.090



Orange filter for the bulb, for light-sensitive patients [01]	C-000.15.359
Fixation cards with holder for dynamic retinoscopy	C-000.15.360

- ∴ **Patented ParaStop®**
- ∴ **Metal controls**
- ∴ **Dustproof**
- ∴ **Streak or Spot Retinoscope by simply changing the bulb**

HEINE Ophthalmic Diagnostic Sets





Set complete with: **BETA 200 Streak Retinoscope**
1 spare bulb
hard case

BETA 200 Retinoscope Sets	2.5V	3.5V (NiMH)	3.5V 
BETA battery handle	C-034.10.118*		
BETA R rechargeable handle for mains socket		C-034.20.376	C-034.27.376
BETA TR rechargeable handle with plug-in transformer		C-034.20.384	C-034.27.384
BETA NT rechargeable handle and NT 300 charger		C-034.20.420	
 BETA L rechargeable handle and NT 300 charger			C-034.29.420



Set complete with: **BETA 200 S Ophthalmoscope, BETA 200 Streak Retinoscope**
one spare bulb each
hard case

BETA Diagnostic Sets	2.5V	3.5V (NiMH)	3.5V 
BETA battery handle	C-262.10.118*		
BETA R rechargeable handle for mains socket		C-262.20.376	C-262.27.376
BETA TR rechargeable handle with plug-in transformer		C-262.20.384	C-262.27.384
BETA NT rechargeable handle and NT 300 charger		C-262.20.420	
 BETA L rechargeable handle and NT 300 charger			C-262.29.420


*Available with large battery handle X-001.99.120 (doesn't fit into hard cases). Please specify with your order.

: Rechargeable handle with charge status display in the bottom insert.

A company makes history and gives a name to outstanding product quality: HEINE.

Founded in 1946, family owned, rich in tradition and today international market leader in a complete range of diagnostic instruments. HEINE products meet all international standards (ISO/CE) and represent the leading edge in precision and ergonomic design. Ongoing dedication to research and development has created a solid and diverse technology base for continued quality and product leadership. Over 500 employees worldwide contribute to this success.

www.heine.com

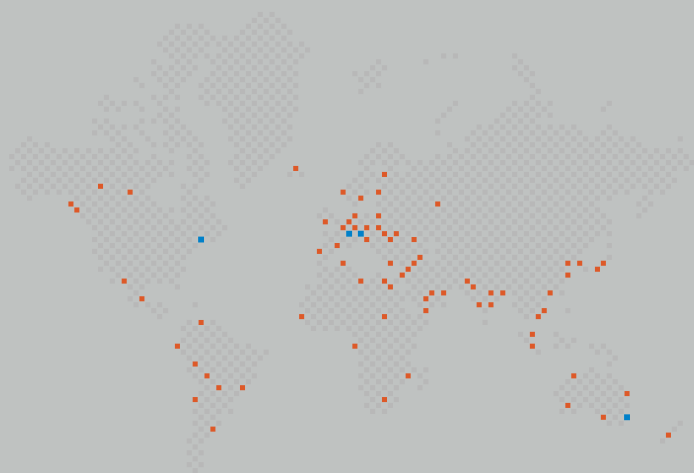

■ HEINE SUBSIDIARY COMPANIES
■ HEINE REPRESENTATIVES

GERMANY

HEINE Optotechnik GmbH & Co. KG
82211 Herrsching
Kientalstr. 7
Tel. +49 (0) 81 52-38 0
Fax +49 (0) 81 52-3 82 02
E-Mail: info@heine.com

NORTH AMERICA

HEINE USA LTD.
10 Innovation Way
Dover, NH 03820
Tel. (603) 7 42-71 03
Fax (603) 7 42-72 17
Toll Free (800) 367-4872
E-Mail: service@heine-na.com

AUSTRALIA

HEINE AUSTRALIA PTY. LTD.
Unit 9, 98 Old Pittwater Road
P.O. Box 717 Brookvale 2100
New South Wales
Tel. +61 (0) 2-99 38 95 00
Fax +61 (0) 2-99 39 23 05
E-Mail: info@heine.com.au

SWITZERLAND

HEINE (Schweiz) AG
Tobeläckerstr. 9
CH-8212 Neuhausen am Rheinfall
Tel. +41 (0) 52-6 72 22 66
Fax +41 (0) 52-6 72 63 77
E-Mail: heine-schweiz@schaffhausen.ch

Your distributor

10/11. A-000.00.100 e

