

# Binocular Loupes

## HIGH MAGNIFICATION AND SHARP, CRISP IMAGES WITH POWERFUL LED ILLUMINATION

HEINE binocular loupes set a benchmark for design, hightech materials and precise manufacturing methods. They provide exceptional image quality and edge-to-edge definition and a large field of view.

Various magnifications (from 2x to 6x), working distances, wearing modes (unique S-Frame, Lightweight Headband and Headband Professional L) and LED LoupeLight illumination systems are available. They provide homogenous, reflex-free illumination even in deep cavities. State-of-the-art LED technology and design!

HEINE Binocular Loupes are used in many disciplines: e.g. General Medicine, Dermatology, Surgery, Cosmetic Surgery, Neurosurgery, ENT, Ophthalmology, Veterinary Medicine, Dentistry etc.



## HEINE Binocular Loupes Introduction

### Definitions:

∴ **Working distance (x)**


The distance, at which the optics of a loupe are sharply-focused.

∴ **Depth of field (y)**

The range over which the image remains sharply-focused.

∴ **Field of view (z)**

The area that can be seen under magnification at the nominal working distance.

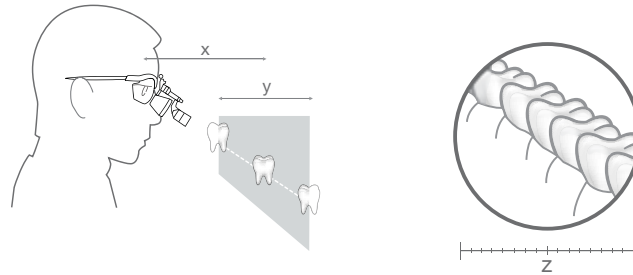

### How to choose the correct Loupe:

∴ **Choose your magnification**

Choose the lowest magnification that meets your needs. The lower the magnification, the larger the field of view and the easier the loupe will be to use.

∴ **Choose your optics**

**HR:** High-Resolution achromatic optics with excellent optical qualities.

**HRP:** High Resolution Prismatic optics with excellent optical qualities for 3.5x magnification and higher.

**HR-C:** Compact and lightweight optics with good optical qualities.

**C:** Economy.

∴ **Choose your working distance**

This is an individual choice depending on the examiner's size and working position. The longer the working distance, the bigger the field of view.

∴ **Choose your mounting option**

Spectacle frame or Headband.

∴ **Add illumination**

Co-axial bright light is especially-important when using magnification.

HEINE Binocular Loupes

	HR	HRP	HR-C	C
				 C2.3K C 2.3
<b>Magnification</b>				
2x	<input type="checkbox"/>			
2.3x				<input type="checkbox"/>
2.5x	<input type="checkbox"/>		<input type="checkbox"/>	
3.5x		<input type="checkbox"/>		
4x		<input type="checkbox"/>		
6x		<input type="checkbox"/>		
<b>Working distance</b>				
250 mm (with close-up lenses)	<input type="checkbox"/>	--	--	--
340 mm	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
420 mm	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	--
450 mm	--	--	--	<input type="checkbox"/>
520 mm	<input type="checkbox"/>	--	--	--
<b>Mounting options</b>				
S-Frame 	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	--
Lightweight Headband 	--	--	<input type="checkbox"/>	<input type="checkbox"/>
Professional Headband 	<input type="checkbox"/>	<input type="checkbox"/>	--	--
Simple Metal Frame 	--	--	--	<input type="checkbox"/>
<b>Accessories</b>				
Close-up Lenses*	<input type="checkbox"/> (for 340mm)	--	--	--
Protective Lenses	<input type="checkbox"/>	<input type="checkbox"/>	--	<input type="checkbox"/>
<b>LED LoupeLight (Pages 129 / 130)</b> 	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	--
<b>Features</b>				
i-View Function**	<input type="checkbox"/>	<input type="checkbox"/>	--	--
Modular	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	--
Water-resistant	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	--
Achromatic Optics	<input type="checkbox"/>	<input type="checkbox"/>	--	--
Optional Prescription Correction	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	--
Page	118-121	122-123	124	125

\* To reduce the working distance to 250 mm.

\*\* i-View function: Provides any angle of view and flips up the optics independently of the illumination beam.



## HEINE HR® 2x High Resolution Binocular Loupes

### with i-View loupe mount



⚡ Multi-coated Achromatic Optics

⚡ Exceptional Focal Range

⚡ Crisp, High-Resolution Image

⚡ Unique i-View Loupe Mount

### HR 2x Optics with i-View

- ⚡ **Achromatic optics** for crisp, bright images with colour-correction.
- ⚡ **Lightweight** (51 g).
- ⚡ **Water-resistant.**
- ⚡ **Independent adjustment** of left and right optics for perfect PD adjustment.
- ⚡ **Largest field of view and greatest depth of field.**
- ⚡ **i-View loupe mount** provides any angle of view and flips up the optics independently of the optional LED LoupeLight.

### HR 2x Specifications

Working distance in mm	340	420	520
Field of view in mm	130	190	250
Depth of field in mm	130	200	290



Optional clip-in correction frame

### Mounting Options:

#### ⚡ HEINE S-Frame®

- ⚡ **Maximum comfort** with a modern, ergonomic design.
- ⚡ **Malleable temples and nose pad** can be customized to fit any face perfectly.
- ⚡ **Innovative materials** – Carbon-Kinetium reduces weight (24 g).
- ⚡ **Tough polycarbonate protective lenses.**
- ⚡ **Optional clip-in correction frame** for prescription lenses.



Professional L Headband with S-Guard

#### ⚡ HEINE Professional L® Headband with S-Guard

- ⚡ **S-Guard splash protection.**
- ⚡ **Ergonomic nose pad** for comfort during long procedures.
- ⚡ **Fully-adjustable for height and width to fit any head.**
- ⚡ **Hygienic and comfortable leather cushions.**
- ⚡ **Optional clip-in correction Frame** for prescription lenses.



Professional L Headband without S-Guard

#### ⚡ HEINE Professional L® Headband without S-Guard

- ⚡ **Ideal for spectacle-wearers in combination with their own spectacles.**
- ⚡ **Ergonomic nose pad** for comfort during long procedures.
- ⚡ **Fully-adjustable for height and width to fit any head.**
- ⚡ **Hygienic and comfortable leather cushions.**

HEINE HR® 2x High Resolution Binocular Loupes

with i-View loupe mount

HR 2x Set with i-View loupe mount



Set A



Set B  
with S-Guard



Set C  
without  
S-Guard

Working distance in mm	340	420	520
<b>Set A:</b> HR 2x optics with i-View on S-Frame Protective lenses 2 sterilisable swivel levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-000.32.375	C-000.32.376	C-000.32.377
<b>Set B:</b> HR 2x optics with i-View on Professional L Headband S-Guard splash protection Protective lenses 2 sterilisable swivel levers Cleaning fluid Cleaning cloth	C-000.32.378	C-000.32.379	C-000.32.380
<b>Set C (without S-Guard):</b> HR 2x optics with i-View on Professional L Headband Protective lenses 2 sterilisable swivel levers	C-000.32.878	C-000.32.879	C-000.32.880

For details of accessories please see pages 127–128 and for LED LoupeLight see pages 129–130.



## HEINE HR® 2.5x High Resolution Binocular Loupes

### with i-View loupe mount



⚡ Multi-coated Achromatic Optics

⚡ Exceptional Focal Range

⚡ Crisp, High-Resolution Image

⚡ Unique i-View Loupe Mount

### HR 2.5x Optics with i-View

- ⚡ **Achromatic optics** for crisp, bright images with colour-correction.
- ⚡ **Lightweight** (51 g).
- ⚡ **Water-resistant.**
- ⚡ **Independent adjustment** of left and right optics for perfect PD adjustment.
- ⚡ **High magnification** with a large field of view and excellent depth of field.
- ⚡ **i-View loupe mount** provides any angle of view and flips up the optics independently of the optional LED LoupeLight.

### HR 2.5x Specifications

Working distance in mm	340	420	520
Field of view	90	130	180
Depth of field	110	180	260



Optional clip-in correction frame

### Mounting Options:

#### ⚡ HEINE S-Frame®

- ⚡ **Maximum comfort** with a modern, ergonomic design.
- ⚡ **Malleable temples and nose pad** can be customized to fit any face perfectly.
- ⚡ **Innovative materials** – Carbon-Kinetium reduces weight (24 g).
- ⚡ **Tough polycarbonate protective lenses.**
- ⚡ **Optional clip-in correction frame** for prescription lenses.



Professional L Headband with S-Guard

#### ⚡ HEINE Professional L® Headband with S-Guard

- ⚡ **S-Guard splash protection.**
- ⚡ **Ergonomic nose pad** for comfort during long procedures.
- ⚡ **Fully-adjustable for height and width to fit any head.**
- ⚡ **Hygienic and comfortable leather cushions.**
- ⚡ **Optional clip-in correction Frame** for prescription lenses.



Professional L Headband without S-Guard

#### ⚡ HEINE Professional L® Headband without S-Guard

- ⚡ **Ideal for spectacle-wearers in combination with their own spectacles.**
- ⚡ **Ergonomic nose pad** for comfort during long procedures.
- ⚡ **Fully-adjustable for height and width to fit any head.**
- ⚡ **Hygienic and comfortable leather cushions.**

HEINE HR® 2.5x High Resolution Binocular Loupes  
with i-View loupe mount

HR 2.5x Sets with i-View loupe mount



Set A



Set B  
with S-Guard



Set C  
without  
S-Guard

Working distance in mm	340	420	520
<b>Set A:</b> HR 2.5x optics with i-View on S-Frame Protective lenses 2 sterilisable swivel levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-000.32.355	C-000.32.356	C-000.32.567
<b>Set B:</b> HR 2.5x optics with i-View on Professional L Headband S-Guard splash protection Protective lenses 2 sterilisable swivel levers Cleaning fluid Cleaning cloth	C-000.32.365	C-000.32.366	C-000.32.367
<b>Set C (without S-Guard):</b> HR 2.5x optics with i-View on Professional L Headband Protective lenses 2 sterilisable swivel levers	C-000.32.865	C-000.32.866	C-000.32.867

For details of accessories please see pages 127–128 and for LED LoupeLight see pages 129–130.

## HEINE HRP® High Resolution Prismatic Binocular Loupes with i-View loupe mount



- ⌘ Multi-coated Achromatic Optics
- ⌘ Prismatic for highest magnification
- ⌘ Crisp, High-Resolution Image
- ⌘ Unique i-View Loupe Mount


### HRP Optics with i-View

- ⌘ **Achromatic optics** for crisp, bright images from edge to edge.
- ⌘ **Lightweight** (85 g).
- ⌘ **Water-resistant.**
- ⌘ **Independent adjustment** of left and right optics for perfect PD adjustment.
- ⌘ **Highest magnification:** Choice of 3.5x, 4x and 6x.
- ⌘ **i-View loupe mount** provides any angle of view and flips up the optics independently of the optional LED LoupeLight.

### HRP Specifications

Working distance in mm	HRP 3.5x/420	HRP 4x/340	HRP 6x/340
Field of view in mm	65	50	40
Depth of field in mm	60	40	30

### Mounting Options:

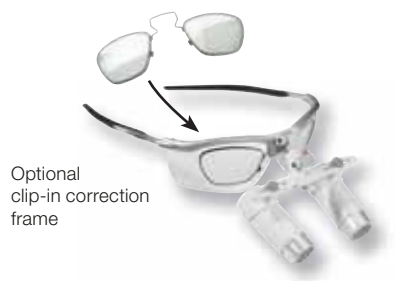

Optional clip-in correction frame

### ⌘ HEINE S-Frame®

- ⌘ **Maximum comfort** with a modern, ergonomic design.
- ⌘ **Malleable temples and nose pad** can be customized to fit any face perfectly.
- ⌘ **Innovative materials** – Carbon-Kinetium reduces weight (24 g).
- ⌘ **Tough polycarbonate protective lenses.**
- ⌘ **Optional clip-in correction frame** for prescription lenses.



Professional L Headband with S-Guard

### ⌘ HEINE Professional L® Headband with S-Guard

- ⌘ **S-Guard splash protection.**
- ⌘ **Ergonomic nose pad** for comfort during long procedures.
- ⌘ **Fully-adjustable for height and width to fit any head.**
- ⌘ **Hygienic and comfortable leather cushions.**
- ⌘ **Optional clip-in correction Frame** for prescription lenses.



Professional L Headband without S-Guard

### ⌘ HEINE Professional L® Headband without S-Guard

- ⌘ **Ideal for spectacle-wearers in combination with their own spectacles.**
- ⌘ **Ergonomic nose pad** for comfort during long procedures.
- ⌘ **Fully-adjustable for height and width to fit any head.**
- ⌘ **Hygienic and comfortable leather cushions.**

HEINE HRP® High Resolution Prismatic Binocular Loupes  
with i-View loupe mount



Set A



Set B  
with S-Guard



Set C  
without  
S-Guard

HRP Sets with i-View

Working distance in mm	HRP 3.5x/420	HRP 4x/340	HRP 6x/340
<b>Set A:</b> HRP optics with i-View on S-Frame Protective lenses 2 sterilisable swivel levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-000.32.430	C-000.32.431	C-000.32.432
<b>Set B:</b> HRP optics with i-View on Professional L Headband S-Guard splash protection Protective lenses 2 sterilisable swivel levers Cleaning fluid Cleaning cloth	C-000.32.440	C-000.32.441	C-000.32.442
<b>Set C (without S-Guard):</b> HRP optics with i-View on Professional L Headband Protective lenses 2 sterilisable swivel levers	C-000.32.840	C-000.32.841	C-000.32.842

For details of accessories please see pages 127–128 and for LED LoupeLight see pages 129–130.



**HEINE HR-C® 2.5x High Resolution Binocular Loupes**  
**with rigid loupe mount**

**HR-C 2.5x Optics**




- ∴ **Multi-Coated optics** for crisp, bright images.
- ∴ **Super-compact lightweight** (17g).
- ∴ **Water-resistant.**
- ∴ **Independent adjustment** of left and right optics for perfect PD adjustment.

- ∴ **Multi-coated Optics**
- ∴ **Compact and lightweight**

**HR-C 2.5x Specifications**

<b>Working distance in mm</b>	<b>340</b>	<b>420</b>
Field of view in mm	<b>60</b>	<b>70</b>
Depth of field in mm	<b>110</b>	<b>140</b>

**Mounting Option:**

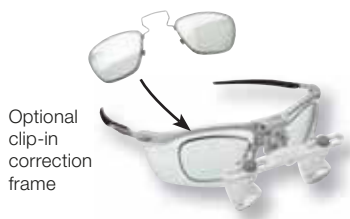

Optional clip-in correction frame

∴ **HEINE S-Frame®**

- ∴ **Maximum comfort** with a modern, ergonomic design.
- ∴ **Malleable temples and nose pad** can be customized to fit any face perfectly.
- ∴ **Innovative materials** – Carbon-Kinetium reduces weight (24g).
- ∴ **Tough polycarbonate protective lenses.**
- ∴ **Optional clip-in correction frame** for prescription lenses.

The HEINE Lightweight Headband is available as an option. See page 127 for ordering information.

**HR-C 2.5x Set with rigid loupe mount**



<b>Working distance in mm</b>	<b>340</b>	<b>420</b>
HR-C 2.5x optics	<b>C-000.32.290</b>	<b>C-000.32.291</b>
S-Frame		
2 sterilisable swivel levers		
Retaining cord		
Cleaning fluid		
Cleaning cloth		
Carrying case		

For details of accessories please see pages 127–128 and for LED LoupeLight see pages 129–130.

## HEINE C 2.3® Binocular Loupes

### The cost-effective all-purpose loupe

Ideal for dental laboratories and general medical examination.



C 2.3

- **Crisp, clear images.** Multi-coated, scratch-resistant silica glass lenses produce crisp, clear images.
- **Flexibility for precise adjustment.** Independently adjustable PD controls ensure precise alignment of the optics.
- **Super lightweight and comfortable – only 46 g.** Flexible coil temples wrap-around to fit the user's anatomy securely and comfortably. Spring-loaded temples will adjust to the user's anatomy. Lightweight metal frame.
- **The ideal all-purpose loupe.** The feather-light weight, low-profile design and low cost of the C 2.3 loupe make it the ideal all-purpose loupe.
- **Low profile design.** Compact, low profile optics allow the user to easily look over the loupes for non-magnified view.
- **Working distance.** 340mm or 450mm (13 or 18 inches).
- Also available as C2.3K on Lightweight Headband.



C 2.3K

Specifications	C2.3x/340mm	C2.3x/450mm
Magnification	2.3x	2.3x
Working distance	340mm	450mm
Field of view dia.	110mm	130mm
Weight	46g	46g
Depth of Field	120mm	190mm

User's refractive error	Effective working distance C2.3x/340mm	Effective working distance C2.3x/450mm
+2	440mm	670mm
+1	390mm	560mm
0	340mm	450mm
-1	300mm	420mm
-2	270mm	370mm

#### Important information for those who wear glasses:

The C 2.3K Series loupe can also be worn without the use of spectacles. Should you prefer to examine without glasses, the charts illustrate how your refractive error will affect the loupe working distance. If you are unsure of your prescription, please consult your eye care provider.

#### C2.3/C 2.3K Binocular Loupes

C2.3/340 Loupe in case	C-000.32.039
C2.3/450 Loupe in case	C-000.32.202
C2.3K/340 Headband Loupe in case	C-000.32.230
C2.3K/450 Headband Loupe in case	C-000.32.231

## Accessories



C2.3/C2.3K Binocular Loupe Accessories	
Protective case for the C 2.3 Loupes and accessories	C-000.32.038
Detachable protective lenses for C 2.3 and C 2.3K. To protect the optics. Can be sterilized by conventional methods.	C-000.32.097



HEINE Binocular Loupes

Ordering Information



 NEW

**Optics with i-View loupe mount for S-FRAME in carrying case**

Working distance in mm	340 mm	420 mm	520 mm
HR 2x	C-000.32.555	C-000.32.556	C-000.32.557
HR 2.5x	C-000.32.595	C-000.32.596	C-000.32.597
HRP 3.5x HRP 4x HRP 6x	C-000.32.721 C-000.32.722	C-000.32.720	

**Optics with i-View and S-Guard for PROFESSIONAL L HEADBAND**

Working distance in mm	340 mm	420 mm	520 mm
HR 2x	C-000.32.558	C-000.32.559	C-000.32.560
HR 2.5x	C-000.32.585	C-000.32.586	C-000.32.587
HRP 3.5x HRP 4x HRP 6x	C-000.32.711 C-000.32.712	C-000.32.710	

**Optics with i-View for PROFESSIONAL L HEADBAND (without S-Guard)**

Working distance in mm	340 mm	420 mm	520 mm
HR 2x	C-000.32.730	C-000.32.731	C-000.32.732
HR 2.5x	C-000.32.733	C-000.32.734	C-000.32.735
HRP 3.5x HRP 4x HRP 6x	C-000.32.737 C-000.32.738	C-000.32.736	

**Optics with rigid loupe mount for S-FRAME and LIGHTWEIGHT HEADBAND**

Working distance in mm	340 mm	420 mm	520 mm
HR-C 2.5x	C-000.32.265	C-000.32.266	

HEINE Binocular Loupes Accessories


Ordering Information



Mounting options	
<b>S-Frame*</b>	
For HR-, HRP-, HR-C Binocular Loupes incl. retaining cord and cleaning fluid	C-000.32.302
* with protective lenses, large (lg).	
<b>Professional L Headband</b>	
For HR and HRP optics	C-000.32.549
<b>Lightweight Headband</b>	
For HR-C optics with rigid loupe mount	C-000.32.516



Accessories for optics	
<b>HR Close-up Lenses in case</b>	
HR 2x/340mm and HR 2.5x/340mm to 250mm. The close-up lenses reduce the working distance of the binocular loupes.	
Close-up Lenses for HR 2x/340mm optics	C-000.32.521
Close-up Lenses for HR 2.5x/340mm optics	C-000.32.523
<b>Protective Lenses</b>	
To protect the optics of the HR and the HRP loupes	
Protective Lenses for HR loupes (5 pairs) [01]	C-000.32.537
Protective Lenses for HRP loupes [02]	C-000.32.027
<b>Sterilisable Swivel Levers</b>	
For HR and HRP loupes with i-View [01]	C-000.32.551
For HR-C loupes with rigid loupe mount [02]	C-000.32.524
<b>Carrying Case</b>	
For Binocular Loupe optics or for all sets with S-Frame	C-000.32.552

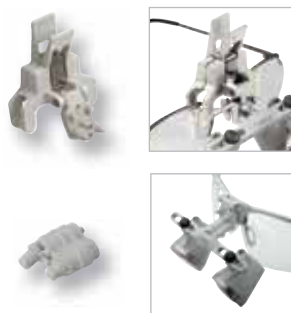

Miscellaneous accessories	
<b>Universal Clip</b>	
For HR-C optics with rigid loupe mount to fit any spectacles	C-000.32.518
<b>Dental Adapter</b>	
For HR-C optics with rigid loupe mount. For a viewing direction from 20° to 35°. It prevents neck and back strain due to poor posture.	C-000.32.528

## HEINE Binocular Loupes Accessories

## Ordering Information

**Accessories for S-Frame**

<b>S-Frame retaining cord</b>	<b>C-000.32.305</b>
-------------------------------	---------------------

**S-Frame replacement part set**

with side rubbers, nose pad (1 piece) and screw

**C-000.32.308****Clip-in correction frame**

for prescription lenses

**C-000.32.309****Protective lenses, large (lg), 1 pair (Standard)****C-000.32.306****Protective lenses, small (sm), 1 pair****C-000.32.307****Cleaning fluid for lenses, pack of 2****C-000.32.542**



## HEINE LED LoupeLight®


### Compact LED LoupeLight, ideal for ENT and dental use




∴ Light

∴ Coaxial shadow-free illumination

∴ 40.000 Lux/25cm distance

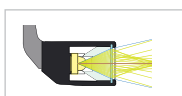

Conventional System [ 01 ]

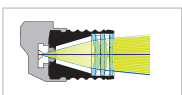

HEINE Precision Lens System [ 02 ]



- ∴ **Advanced LED Technology.** White, homogenous light, ideal for all surgical and dental applications.
- ∴ **Coaxial Illumination.** Shadow-free illumination of cavities, ideal for ENT and dental applications. Perfect co-axial viewing with all compatible HEINE Loupes.
- ∴ **Compact, lightweight (24 g).** Ensures perfect comfort.
- ∴ **Aluminium housing.** Compact and long lasting with advanced heat management. Reduces heat build-up.
- ∴ **Waterproof design.** Easy to clean.
- ∴ **Individual adjustment of the angle of view.** The right angle for every situation.
- ∴ **Optional yellow filter** reduces blue light. Avoids premature hardening of composites.

### Technical Specification

Illumination	<b>40,000 Lux</b> at a working distance of 250mm
Colour Temperature	<b>6,250 Kelvin</b>
LED Operating Life	1W Super-High Output LED rated at <b>50,000 hours operating life</b>
Spot Size	The illuminated field is designed to complement the HR/HRP loupes. <b>81 mm dia. at 420 mm, 95 mm dia. at 500 mm working distance</b>
Angle of Declination	Individually-adjustable for every examination
Operating Time	17 hours (with mPack LL)

- ∴ **Precision multiple lens optical system** consisting of four high-quality lenses. Completely homogenous illumination of the whole field of view and a clearly defined light spot.


## HEINE mPack LL®

### Power Source




- ∴ **Portable power pack.** Does not tie the user to desk, wall or static power supplies.
- ∴ **Li-ion Technology.** Triple capacity compared with current NiCd systems. No memory effect, autonomy for 17 hours.
- ∴ **Fast Charge.** Charges to full capacity in 2 hours.
- ∴ **Charge Status Display.** Always lets you know how much charge you have left.
- ∴ **Stepless light intensity control.** Optimum brightness setting prevents reflexes.
- ∴ **Flexible Charging Options:** EN 50 desk and wall charger or plug-in transformer.
- ∴ **Intelligent.** Can be used as mains power supply even while in "charge" mode.
- ∴ **100-240 V Compatible Transformer** incl. 4 plugs for world wide compatibility.

For the optional EN50 table/wall unit please see page 153.

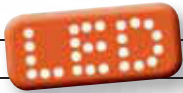

HEINE LED LoupeLight® and Accessories

Ordering Information

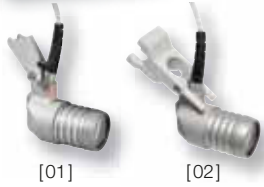


[01]

[02]



[03]

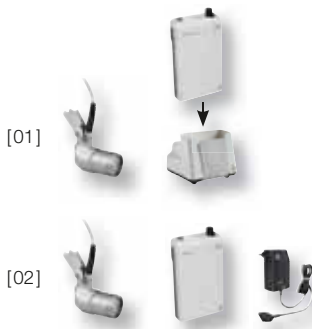

[01]

[02]

LED LoupeLight

LED LoupeLight for HR/HRP Loupes with i-View loupe mount [02]	C-008.32.236
LED LoupeLight for HR-C Loupes with rigid loupe mount [01]	C-008.32.238

Detachable Yellow Filter [03]	C-000.32.241
-------------------------------	--------------

LED LoupeLight with Power Source

	For i-View loupe mount (HR / HRP)	For rigid loupe mount (HR-C)
LED LoupeLight with mPack LL and EN 50 [01]	C-008.32.381	C-008.32.382
LED LoupeLight with mPack LL and plug-in transformer [02]	C-008.32.384	C-008.32.385

LED LoupeLight Sets with HR Binocular Loupes



Magnification and working distance in mm	HR 2.5x/340	HR 2.5x/420	HR 2.5x/520
HR 2.5x optics with i-View on S-Frame LED LoupeLight Detachable Yellow Filter mPack LL Charger EN 50 Protective lenses 2 sterilisable levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-008.32.355	C-008.32.356	C-008.32.357




LED LoupeLight Sets with HRP Binocular Loupes

Magnification and working distance in mm	HRP 3.5x/420	HRP 4x/340	HRP 6x/340
HRP optics with i-View on S-Frame LED LoupeLight Detachable Yellow Filter mPack LL Charger EN 50 Protective lenses 2 sterilisable levers Retaining cord Cleaning fluid Cleaning cloth Carrying case	C-008.32.430	C-008.32.431	C-008.32.432

**A company makes history and gives a name to outstanding product quality: HEINE.**

Founded in 1946, family owned, rich in tradition and today international market leader in a complete range of diagnostic instruments. HEINE products meet all international standards (ISO/CE) and represent the leading edge in precision and ergonomic design. Ongoing dedication to research and development has created a solid and diverse technology base for continued quality and product leadership. Over 500 employees worldwide contribute to this success.

[www.heine.com](http://www.heine.com)

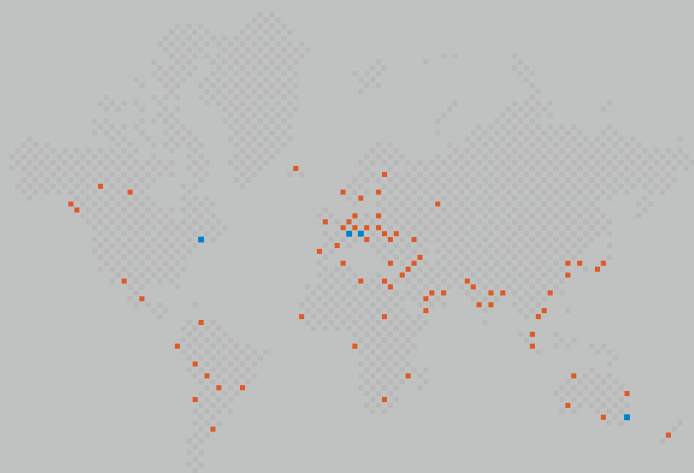

■ HEINE SUBSIDIARY COMPANIES  
■ HEINE REPRESENTATIVES

**GERMANY**

HEINE Optotechnik GmbH & Co. KG  
82211 Herrsching  
Kientalstr. 7  
Tel. +49 (0) 81 52-38 0  
Fax +49 (0) 81 52-3 82 02  
E-Mail: [info@heine.com](mailto:info@heine.com)

**NORTH AMERICA**

HEINE USA LTD.  
10 Innovation Way  
Dover, NH 03820  
Tel. (603) 7 42-71 03  
Fax (603) 7 42-72 17  
Toll Free (800) 367-4872  
E-Mail: [service@heine-na.com](mailto:service@heine-na.com)

**AUSTRALIA**

HEINE AUSTRALIA PTY. LTD.  
Unit 9, 98 Old Pittwater Road  
P.O. Box 717 Brookvale 2100  
New South Wales  
Tel. +61 (0) 2-99 38 95 00  
Fax +61 (0) 2-99 39 23 05  
E-Mail: [info@heine.com.au](mailto:info@heine.com.au)

**SWITZERLAND**

HEINE (Schweiz) AG  
Tobelackerstr. 9  
CH-8212 Neuhausen am Rheinfall  
Tel. +41 (0) 52-6 72 22 66  
Fax +41 (0) 52-6 72 63 77  
E-Mail: [heine-schweiz@schaffhausen.ch](mailto:heine-schweiz@schaffhausen.ch)

Your distributor

10/11. A-000.00.100 e

