

∴ HEINE® iC 2 Funduscope

If the patient cannot come to the instrument, the instrument can simply go to the patient.



LED HQ

LED NOW IN HEINE QUALITY.



Mobile, digital ophthalmoscopy for children and bedridden patients. And everyone else.



✓ **Always with you wherever you are.**
With a weight of just 300g, practical dimensions and cable-free operation, it's ideal for examining children as well as bedridden, injured and disabled patients.

✓ **Comfortable for both the doctor and patient.**
In addition to the conveniently longer working distance for both sides, full adjustment of LED brightness is particularly pleasant for the patient. Needless to say, it corresponds to the highest international standard for photobiological safety.

✓ **Continuous diopter compensation.**
Adjustable from +15 to -15 dioptres, complemented by autofocus with +/-3 dioptres.

NEW!
✓ **Set it down and charge in one step.**
The new iC2 Charger is the world's first inductive charging station to receive approval for a medical product. Say goodbye to charging cables!

✓ **High-quality fundus images:**
Optical precision lens system with 4 aspherical high-performance lenses. Convenient documentation of completed examinations. Excellent imaging quality at up to 12 MP**.

✓ **Field of view up to 34°.**
The angle of view is therefore approximately four times as big as that of conventional direct ophthalmoscopes with 9°.

✓ **Straightforward scanning of the periphery,**
simply by tilting the iC2 back and forth.



✓ **See for yourself what the iC2 app is capable of:**
intuitive localisation of images, continuous magnification, storage of additional patient data, documentation of the current situation and, most notably, provision of images via email.



✓ **A film says more than a thousand words.**
A brief film shows you just how practical, simple and fast fundus photography can be. Scan this QR code or go to heine.com/video/heine_ic2_en_gb/



Examples of fundus images with the HEINE iC2 Funduscope

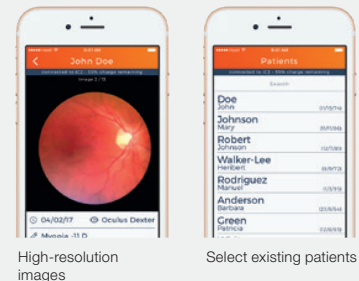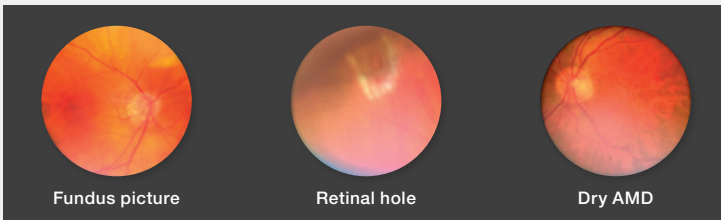

High-resolution images

Select existing patients

Ordering information:

iC2 Charger Set/7 for iPhone 7 and 8

iC2 with mounting case smartphone iC 2/7 for Apple® iPhone® 7 and 8*, iC2 Charger, USB cord with E4-USB medical approved plug-in power supply and case **C-152.40.211**

iC2 Charger Set/6 for iPhone 6/6s

See above, but with mounting case smartphone iC 2/6 for Apple iPhone 6/6s* **C-151.40.211**

iC2 Charger Set/5 for iPhone 5s and SE

See above, but with mounting case smartphone iC 2/5 for Apple iPhone 5s and SE* **C-150.40.211**



▶ Due to local regulations, some products may not be available in your country. For additional product configurations and further HEINE Quality Products please visit heine.com

* The Apple iPhone is not included in the scope of delivery. Apple, the Apple logo and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

** With iPhone from Version 6s

HEINE Optotechnik GmbH & Co. KG, Kientalstr. 7, 82211 Herrsching, Germany
Tel. +49 (0) 8152/38-0, Fax +49 (0) 8152/38-202, E-Mail: info@heine.com, www.heine.com



We reserve the right to change specification without notice.
© 2018 HEINE Optotechnik GmbH & Co. KG. All rights reserved. A-000.00.219 e 10/18.