

Welcome to the  
new standard in  
Digital Dermatoscopy.



Where we combine proven  
HEINE Quality with the  
digital future.




# Simple, flexible and at a fair price.

Skin screening, documentation, data protection – the new HEINE Cube is the missing link in modern dermatoscopy and allows you to make more out of your day.



An innovative  
documentation system  
celebrates its  
world premiere.

THE HEINE CUBE

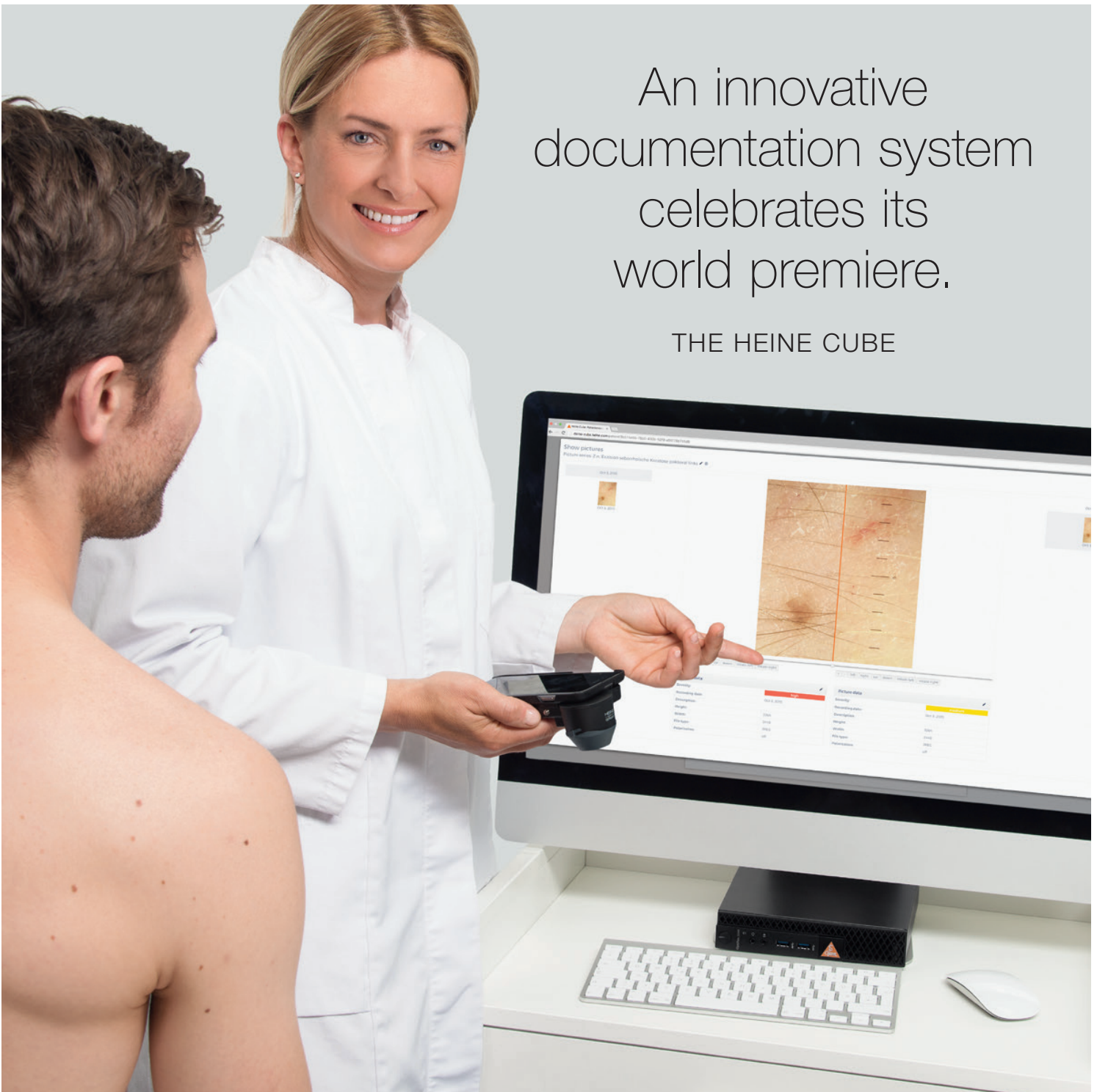

## Simple Set Up.

It's as simple as connecting your new mobile phone or new computer to your wireless network. Unpack the HEINE Cube, switch it on and use the provided access code to connect to the HEINE Cube integrated Wireless Network Hotspot from your PC, Laptop or Tablet web browser. Download the HEINE Cube App to your mobile phone and you are ready to go.

## Fast Image Capture.

With the iPhone\* in combination with the HEINE iC1 you can now capture high quality clinical or dermatological images fast and cable-free. You can then easily assign the image to a body map, make notes and automatically update the patient file. The high quality image (up to 12 megapixels\*\*) is then transmitted wirelessly and in real time to the HEINE Cube. The image is immediately available for you or other users on their connected devices.

## Easy to use.


You can now conveniently work with the information from your choice of device such as a Tablet, Laptop, or your office PC through the Web Browser Interface. Easily compare images with 40x true magnification, add notes and further information, or print out the Patient Information as a PDF. The HEINE Cube integrated Web Browser Interface is platform-independent, so it works on any operating system such as Mac (Apple) or PC (Windows). GDT\*\*\* and DICOM\*\*\*\* are used as standard data transfer interfaces.




### HEINE Cube

Integrated Wireless Network Hotspot with WPA 2 encryption for stand-alone use and / or as network interface to integrate into your practice network.

**HEINE iC1**  
with HEINE Precision Optics combined with the high resolution iPhone\* camera for the best dermatoscopic images.

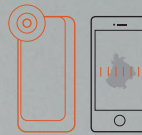

**Mounting cases smartphone** are available for iPhone\* from version 5s upwards or iPod\* touch 6th generation. The iPhone/iPod\* is not included in the scope of delivery.

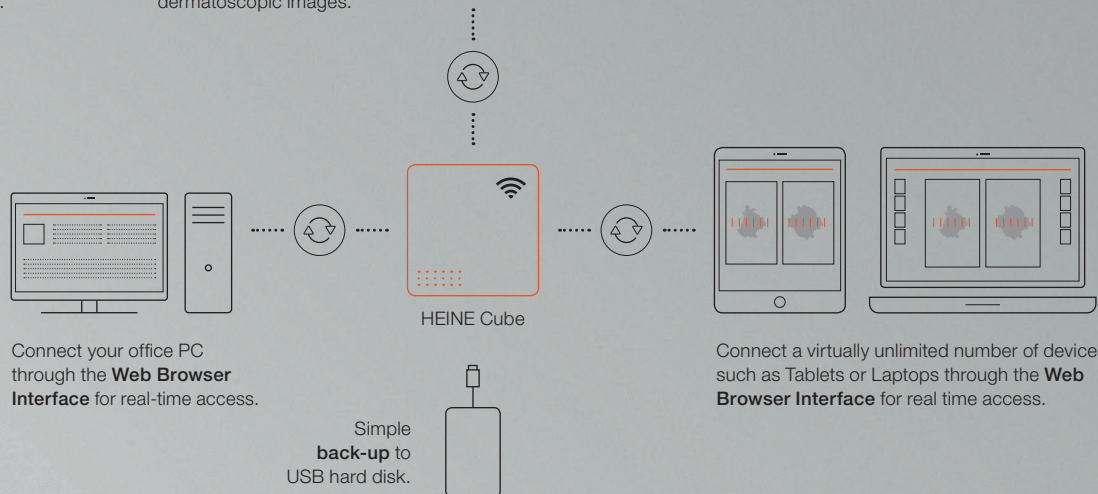

\* Apple, the Apple logo, iPhone and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. \*\* with iPhone from Version 6s  
\*\*\* GDT [German Standard for Device Data Transfer enabling system-independent data transfer between medical devices and external programs such as Practice Management of Electronic Health Record Systems] \*\*\*\* DICOM [Digital Imaging and Communications in Medicine]



For more information about  
the HEINE Cube, please visit:  
[www.heine.com/en/digital](http://www.heine.com/en/digital)

Or call us: +49(0)81 52/3 85-00



HEINE Optotechnik GmbH & Co. KG  
Kientalstr. 7, 82211 Herrsching, Germany  
[www.heine.com](http://www.heine.com)



We reserve the right to change specification without notice.  
© 2018 HEINE Optotechnik GmbH & Co. KG. All rights reserved. A-000.00.270 e 09/18.